Anti-Carbonic Anhydrase IV / Car4 / CA4 antibody

Catalog Number: 100708



Product name

Anti-Carbonic Anhydrase IV / Car4 / CA4 antibody

Immunogen

<u>Human Carbonic Anhydrase IV / Car4 / CA4 (His Tag) recombinant protein</u>

Specificity

Human Carbonic Anhydrase IV / CA4

No cross-reactivity in ELISA with Human CA2;
Human CA3; Human CA5A; Human CA5B; Human
CA8; Human CA9; Human CA10; Human CA12;
Human CA13; Human CA14

Antibody description

Rabbit monoclonal to Carbonic Anhydrase IV / Car4 / CA4

Preparation

This antibody was obtained from a rabbit immunized with purified, recombinant Human Carbonic Anhydrase IV / CA4 (rh Carbonic Anhydrase IV / CA4; NP_000708.1; Met1-Lys283).

Formulation

0.2 µm filtered solution in PBS with 5% trehalose

Storage

This antibody can be stored at 2°C-8°C for one month without detectable loss of activity. Antibody products are stable for twelve months from date of receipt when stored at -20°C to -80°C .

Preservative-Free.

Sodium azide is recommended to avoid contamination (final concentration 0.05%-0.1%). It is toxic to cells and should be disposed of properly. Avoid repeated freeze-thaw cycles.

Clonality

Monoclonal

Ig Type

Rabbit IgG

Applications

ELISA, IHC-P


Dilutions

ELISA: $0.1-0.2 \mu g/mL$

This antibody can be used at 0.1-0.2 μ g/mL with the appropriate secondary reagents to detect Human Carbonic Anhydrase 4 / CA4. The detection limit for Human Carbonic Anhydrase 4 / CA4 is approximately 0.039 ng/well.

IHC-P: 1-10 μg/mL

Validations

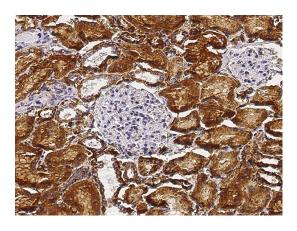


Carbonic Anhydrase IV / CA4 Antibody, Rabbit MAb, Immunohistochemistry

Immunochemical staining of human CA4 in human kidney with rabbit monoclonal antibody (5 μ g/mL, formalin-fixed paraffin embedded sections). Positive staining was localized to renal tubules.

Anti-Carbonic Anhydrase IV / Car4 / CA4 antibody



Catalog Number: 100708

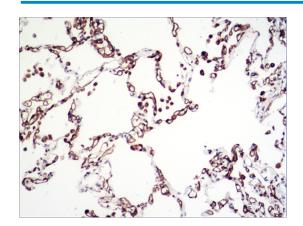

Carbonic Anhydrase IV / CA4 Antibody, Rabbit MAb, Immunohistochemistry

Immunochemical staining of human CA4 in human lung with rabbit monoclonal antibody (5 μ g/mL, formalin-fixed paraffin embedded sections). Positive staining was localized to pulmonary capillaries.